


## Tetracycline USP Assay by HPLC- AppNote

## **USP Method for Tetracycline**

This USP Tetracycline Assay is performed and a five injection overlay shows run to run consistency and Precision. The Peak Efficiency is very good and Peak Tailing guidelines are easily met with RSD values less than 0.4%. This demonstrates a good Alternate Column for your USP Method.

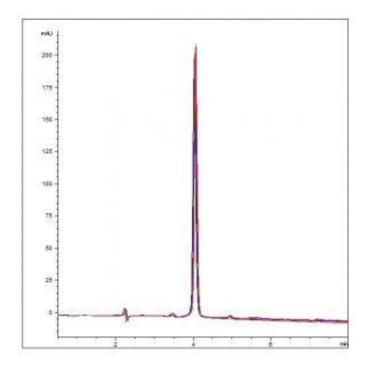

## Peak:

Tetracycline

## **Method Conditions:**

Column: Cogent RP C18 PE™, 5µm, 100Å

**Catalog No.:** <u>84518-15P</u> **Dimensions:** 4.6 x 150mm

**Mobile Phase:** 

A: DI Water 0.1% Phosphoric Acid

B: Acetonitrile

Time (minutes)

0	15
7.5	15
7.6	40
10	15

Injection vol.: 10µL

Flow rate: 1.0 mL / minute Detection: UV @ 280nm

Sample Preparation: Standard solution: 100µg / mL of USP Tetracycline in Mobile Phase A

Most recently appeared in Pharmacopeial Forum: Volume No. 43(6) Page Information:

USP43-NF38 – 1827 USP42-NF37 – 1787 USP41-NF36 – 1690

Note: Tetracycline is a broad spectrum antibiotic produced by the Streptomyces genus of Actinobacteria. It exerts a bacteriostatic effect on bacteria by binding reversible to the bacterial 30S ribosomal subunit and blocking incoming aminoacyl tRNA from binding to the ribosome acceptor site. It also binds to some extent to the bacterial 50S ribosomal subunit and may alter the cytoplasmic membrane causing intracellular components to leak from bacterial cells.



Printed from the Chrom Resource Center
Copyright 2025, All Rights Apply
MicroSolv Technology Corporation
9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900
Fax: (910) 769-9435
Email: customers@mtc-usa.com
Website: www.mtc-usa.com