
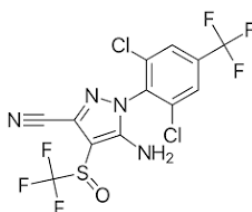


Fipronil Analyzed with HPLC - AppNote

A Reproducible Method for a Phenylpyrazole Insecticide

A rapid, sensitive, and reproducible method has been developed for this Insecticide. The data below, *(an overlay of 10 chromatograms)* illustrates how the compound can be adequately retained and detected using this straightforward method.

A phenyl ring in the column stationary phase provides strategic use of π - π Interaction with the analyte making possible the use of a very simple, Mass Spec-friendly mobile phase with formic acid as an additive.

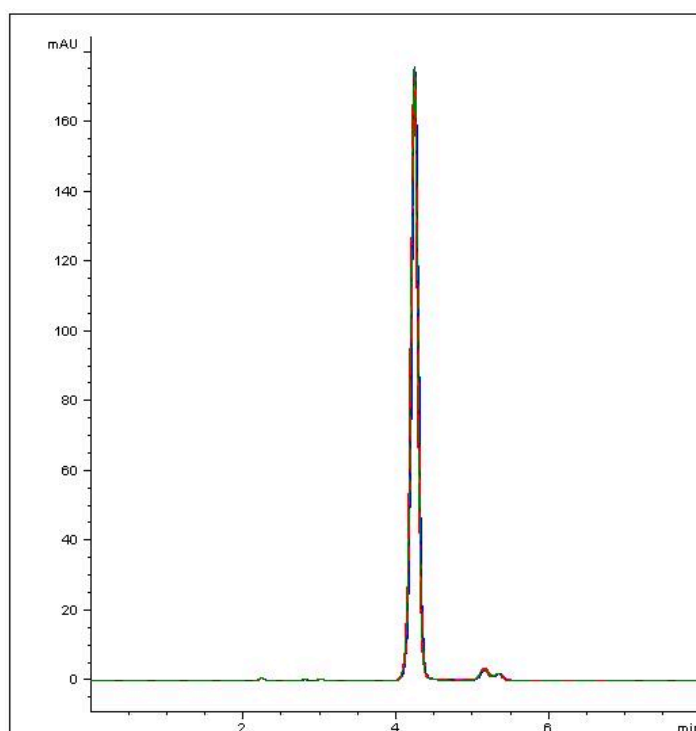

Peak:

Fipronil

Method Conditions

Column: Cogent Phenyl Hydride™, 4 μ m, 100Å

Catalog No.: [69020-15P](#)

Dimensions: 4.6mm x 150mm

Mobile Phase: 60:40 Acetonitrile / DI Water with 0.1% Formic Acid

Injection vol.: 1µL

Flow rate: 1.0mL / minute

Detection: UV @ 254nm

Sample Preparation: 0.01% Solution in 50:50 Acetonitrile / DI Water

%RSD: <0.1%

t₀: 1.5 Minutes

K': 2.5

Notes: Fipronil is a Phenylpyrazole insecticide and Acaricide. It was originally used in agriculture as seed coatings, but is commonly used as flea and tick control in dogs and cats. Fipronil is believed to act as a noncompetitive blocker of GABA-gated chloride channels.

Note 2: Capacity is determined using the following equation: $k = (t_R - t_0) / t_0$

- t_R = Retention Time of an Analyte Peak
- t_0 = Retention Time of non-Retained Peak



Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com